

Pelvic Stabilisation Device


TPB – Trauma Pelvic Binder




**100%
Radiolucent**




**Uniform
compression**

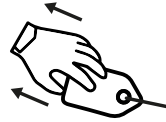

**Easy one-handed
closure**

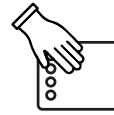

**Can be handled by
a single professional**

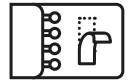

**Removable X-ray
detectable strap**

TPB Trauma Pelvic Binder enables circumferential, uniform and controlled compression over the pelvic region quickly and easily.

Designed for patients with suspected pelvic injury, TPB allows simultaneous circumferential closure that compresses the pelvis precisely, promoting early stabilisation and reducing the risk of haemorrhage. Its design prevents overpressure, ensuring safe and effective application.

Thanks to its pulley system, the device can be easily positioned and adjusted with one hand, allowing it to be used by a single trained healthcare professional, always in accordance with established local protocols.

TPB is 100% radiolucent and compatible with X-rays, magnetic resonance imaging (MRI) and computed tomography (CT). It also incorporates a removable, X-ray detectable strip that facilitates identification of the device in the patient and improves the interpretation of radiological images.

The TPB belt combines innovation, ease of use and precision, making it an essential tool for the management of patients with pelvic injuries. It complies with the *American College of Surgeons'* best practice guidelines for the treatment of pelvic fractures with associated haemorrhage.

Pelvic Stabilisation Device

TPB – Trauma Pelvic Binder

Radiological compatibility

Compatible with X-ray, magnetic resonance imaging (MRI) and computed tomography (CT) procedures

Integrated pulley system

Distributes closure force evenly, combining mechanical precision and ease of use

Ergonomic

Designed for faster and more precise placement



Compact

Compact design to optimise space and portability

Application History Label

Label for easily noting the date and time of application

Product specifications

| Reference | Product description | Dimensions | Color |
|----------------------|----------------------------------|---------------------------------------|----------|
| EME-100-001 (Orange) | Pelvic Stabilisation Device TPB | 140 cm x 20 cm | ● Orange |
| EME-100-002 (Black) | Pelvic Stabilisation Device TPB | 140 cm x 20 cm | ● Black |
| Packaging (1 unit) | Dimensions: 30 cm x 20 cm x 5 cm | MOQ (Minimum Order Quantity) 10 units | |

Instructions for use



© 2026 All rights reserved | IES MEDICAL® (International Emergency Services, S.L.)



CAT_TP_B_EN_01